

Universal
Specimen
Preservation

eSwab

Molecular &
Culture Transport

Collect • Snap • Elute



1 ml Liquid Amies
with Patented
Flocked Swab*


Sample rapidly elutes
into liquid phase

Provides 1ml for multitest analysis

Preserves Nucleic Acids and
Antigens of bacteria, viruses
and Chlamydia

Maintains viability of aerobic,
anaerobic and
fastidious bacteria

Certified
DNase
and
RNase free

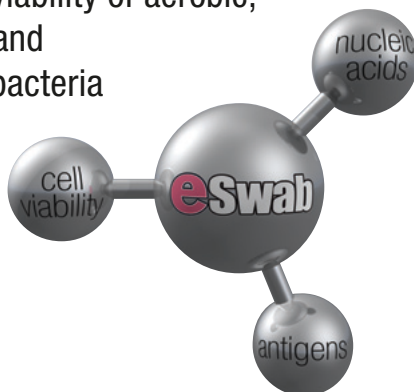

Collect



Snap




Elute



Universal open platform for bacterial culture and the preservation of DNA, RNA and antigens of diverse range of microorganisms. Features unique molded breakpoint with colored indication mark and swab capture cap.

www.copaninnovation.com

* Patented Worldwide: WO 2004/086979

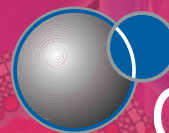

COPAN
innovation




eSwab molecular & culture transport


Catalog No.	Copan eSwab Collection Kit - Product Descriptions	Pack Size	Color Code	Sampling Sites	Capture Cap Feature
480CE	Sterile single use sample collection pack containing: Polypropylene screw-cap tube with internal conical shape filled with 1ml of Liquid Amies Medium. One regular size applicator swab with flocked nylon fiber tip.	50 kits per shelf pack 10 x 50 kits per box	Pink	Nose, throat, vagina, rectum, faces and wounds	YES
481CE	Sterile single use sample collection pack containing: Polypropylene screw-cap tube with internal conical shape filled with 1ml of Liquid Amies Medium. One minitip applicator swab with flocked nylon fiber tip.	50 kits per shelf pack 10 x 50 kits per box	Orange	Eye, ear, nasal passages, throat, urogenital tracts	NO
482CE	Sterile single use sample collection pack containing: Polypropylene screw-cap tube with internal conical shape filled with 1ml of Liquid Amies Medium. One pernasal applicator swab with flocked nylon fiber tip.	50 kits per shelf pack 10 x 50 kits per box	Blue	Nasopharynx and pediatric sample collection	NO
483CE	Sterile single use sample collection pack containing: Polypropylene screw-cap tube with internal conical shape filled with 1ml of Liquid Amies Medium. One urethral applicator swab with flocked nylon fiber tip.	50 kits per shelf pack 10 x 50 kits per box	Orange	Urogenital tract	NO
484CE	Sterile single use sample collection pack containing: Polypropylene screw-cap tube with internal conical shape filled with 1ml of Liquid Amies Medium. One pediatric applicator swab with flocked nylon fiber tip.	50 kits per shelf pack 10 x 50 kits per box	Blue	Pediatric sample collection	NO

Note: Other eSwab product codes are available for the use with automated streaking systems for culture

How to Use the eSwab

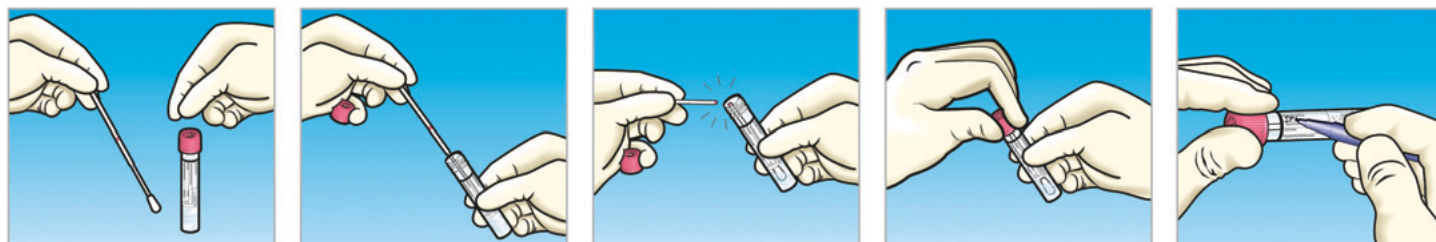

How the Cap Captures the Broken Swab Applicator

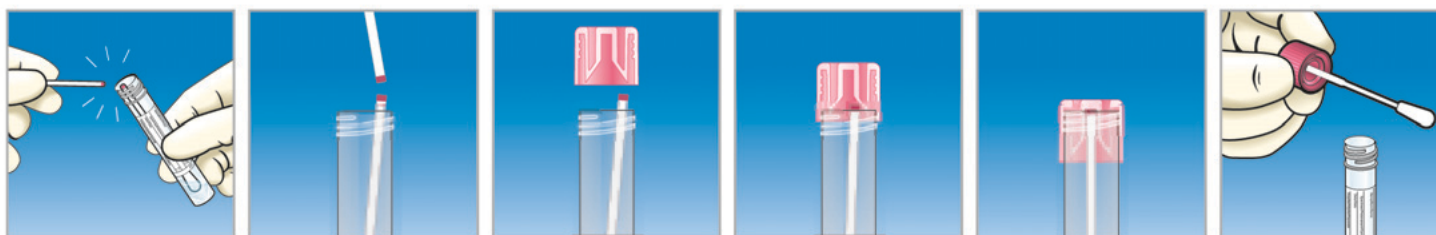


In the testing laboratory, the eSwab cap is unscrewed and removed, the swab applicator stick is securely attached to the cap. This feature allows the operator to conveniently remove the swab and perform various microbiology analyses using the tube cap as a handle to hold and manipulate the swab. Due to the flexibility of the shaft of 481CE, 482CE, 483CE and 484CE the capture cap feature is not applicable, as the broken applicator may not firmly fit into the cap.

K65 01.09

Collection Swab Breakpoint

The swab applicator shaft is broken off at the colored breakpoint indication line into the eSwab tube of transport medium.

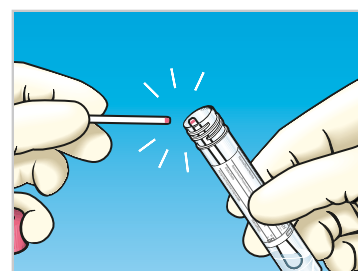

Copan Italia S.p.a. Via F. Perotti, 10, Brescia, 25125, Italy

Tel: +39 030 268 7211 • Fax: +39 030 268 7250

E-mail: info@copanitalia.com Check out our website: www.copaninnovation.com

Copan Flocked Swabs:
Worldwide Patent No.
WO 2004/086979

CE 0344